




**HPQHL-30 Strong Anion Exchange  
Chromatography Resin  
Product Manual**

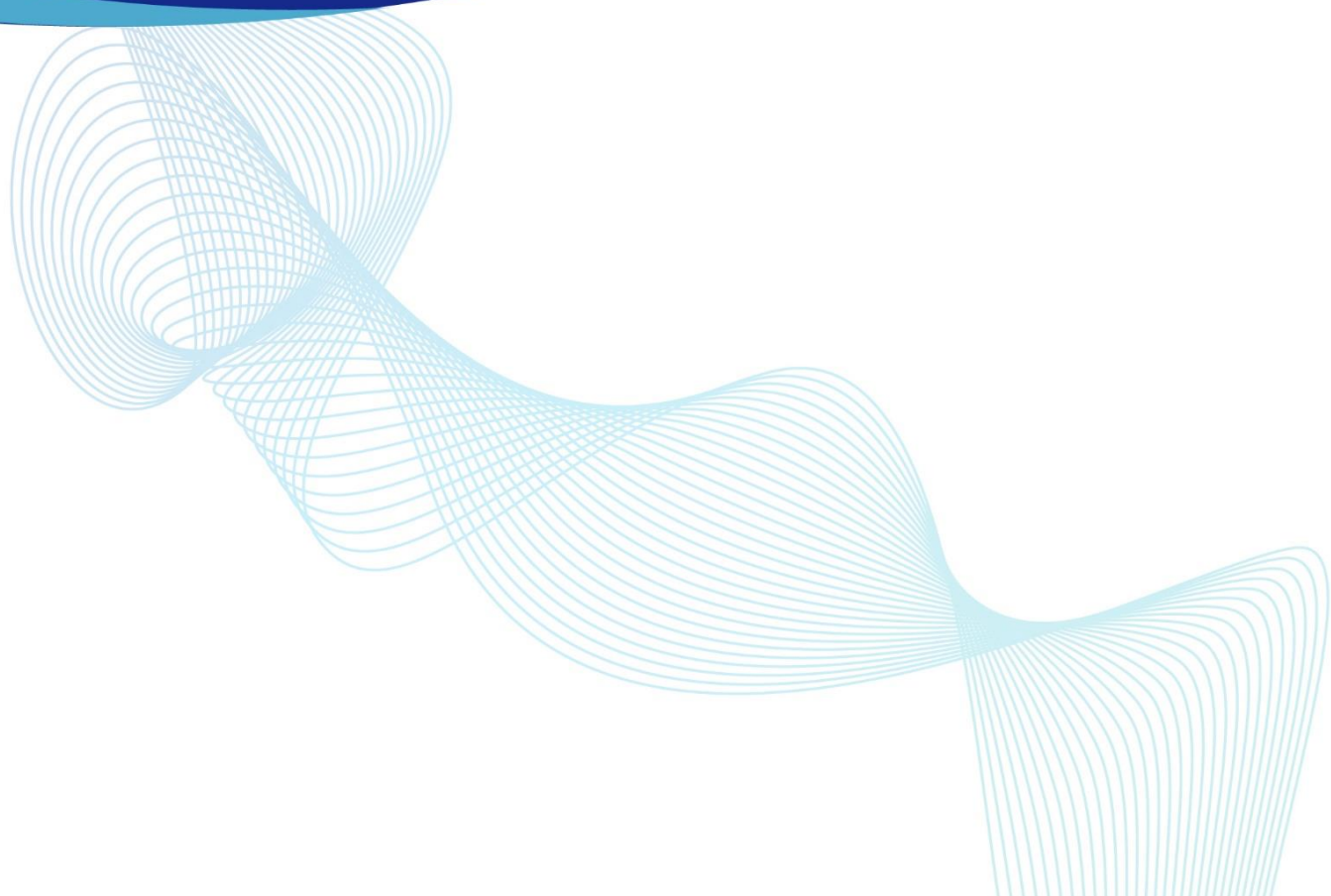

## 1. Product Introduction

Ion exchange chromatography is the most commonly used method for the separation and purification of biomacromolecules. HPQHL-30 strong anion exchange chromatography resin developed by Bogen is based on polymethacrylate microspheres with rigid structure as matrix. The hydrophilicity of the polymer matrix is improved through unique surface modification technology, and the adsorption performance of functional groups is improved with Bogen's unique high-capacity ligand synthesis technology. This resin features the advantages of high resolution, high binding capacity, low non-specific adsorption and excellent biocompatibility, is suitable for the polishing of protein, enzyme, polysaccharide and nuclear acid.

## 2. Product Properties

Parameter	Technical Specification
Particle size range	40±10 μm
Pore size	30 nm
Matrix beads	Polymethacrylate
Function group	Quaternary amino group (-CH <sub>2</sub> N <sup>+</sup> (CH <sub>3</sub> ) <sub>3</sub> )
Function group type	Strong base
Dynamic binding capacity	≥85mg BSA/ml wet gel
Pressure upper limit	1 MPa
pH stability	Resistant to 0.5M NaOH solution cleaning
Storage	4-30 °C (20% ethanol)

## 3. Operation Steps


HPQHL-30 strong anion exchange resin is widely used for polishing of proteins, enzymes, polysaccharides, and nucleic acids. Chromatography operations usually include steps such as equilibration, loading, washing, elution, and regeneration.

Detailed operation steps are as follows:

**Sanitization:** Rinse the column with 0.5 M NaOH at a low flow rate for 30 minutes.

**Equilibration:** Use 10 - 20 CVs equilibration buffer (Buffer A, such as 20 mM Tris-HCl, pH 8.0, the actual buffer system used should be screened and optimized based on the stability and isoelectric point of the target protein) at 250-500 cm/h flow rate to equilibrate the ion exchange column until the conductivity and pH of the effluent remain stable (consistent with the equilibration solution).

**Loading:** The buffer of sample should be as consistent as possible with the equilibration solution. Solid samples can be prepared by dissolving in the equilibration solution; low-concentration sample solutions can be dialyzed with the equilibration solution or add a corresponding amount of salt;



---

high-concentration sample solutions can be diluted with the equilibration solution. To avoid clogging the chromatography column, the sample solution should be centrifuged or micro-filtrated (0.45  $\mu\text{m}$ ). The loading amount is calculated based on the loading capacity of the resin and the content of the target protein in the feed.

**Washing:** After loading the sample, continue to wash with equilibration buffer until the UV value drops to the baseline.

**Elution:** Elute with elution buffer (such as 20 mM PB + 1 M NaCl, pH 8.0). The elution method can be pH gradient elution, linear gradient elution or step gradient elution, and collect the effluent.

**Regeneration:** The column can be washed with 1-2 M NaCl after each campaign to remove strongly bound proteins.

**Cleaning in place:** After the resin has been used 5-10 times (the actual number of uses is related to the type and source of feed materials and process requirements), the resin needs to be cleaned in place:

- (1) For proteins that are strongly bound by ionic bonds, clean with 3-4 CVs of 2 M NaCl;
- (2) For precipitated proteins, hydrophobically bound proteins, and lipoproteins, clean with 3-4 CVs of 0.5 M NaOH;
- (3) For strongly hydrophobically bound proteins, lipoproteins and lipids, clean with 3-4 CVs of 70% ethanol or 30% isopropyl alcohol at a flow rate not higher than the upper pressure limit (when using high-concentration organic solvents, in order to avoid bubbles, a method of gradually increasing the concentration of organic solvent should be used).

**Storage:** Store in 20% ethanol at 4-30°C; the resin in the chromatography column can be washed with 20% ethanol and stored at 4-30°C.

**Other precautions:** During operating and storing the column, avoid the column from drying out or being loosely sealed to prevent air bubbles from entering.



## Duoning Biotechnology Group

■ [marketing@duoningbio.com](mailto:marketing@duoningbio.com)

● [www.duoningbio.com/en](http://www.duoningbio.com/en)

📍 6/F manulife place 348kwun tong road kl, HongKong, PRC.

